

# Pipetman® P

## The Pipetting Standard

Designed and manufactured to provide robust, accurate and precise liquid sample measurements, eight models cover the full range of volumes from 0.2 µL to 10 mL.

### ◎ The standard for robustness

- Built to last. Constructed of stainless-steel and PVDF (a plastic with high resistance to chemical attack).

### ◎ The standard for accuracy and precision


- Unique “dry seal” technology avoids using grease to obtain adequate air-tightness and guarantees long-term performance.
- Thick body (handle) in PVDF minimizes the hand-warming effect during pipetting.
- Direct volume reading in decimal digits helps prevent mistakes in volume setting.



Digital volume indicator

### ◎ Expanded flexibility

- Models P2 and P10 incorporate dual position tip-ejectors. Allows use of Diamond® D10 (short) and DL10 (long) tips – a plus for protection and precision when using microtubes.

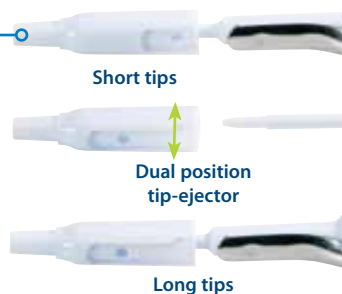

Short tips

Dual position tip-ejector

Long tips

Easy maintenance and cleaning



Pipetman Pipettes				
Model	Diamond Tips	Part Number	Volume (µL)	Price
P2 0.2 - 2 µL	D10 DL10	F144801	Min. 0.2	218,40 €
			0.5	
			Max. 2	
P10 1 - 10 µL	D10 DL10	F144802	Min. 1	218,40 €
			5	
			Max. 10	
P20 2 - 20 µL	D200	F123600	Min. 2	181,20 €
			5	
			Max. 20	
P100 20 - 100 µL	D200	F123615	Min. 20	181,20 €
			50	
			Max. 100	
P200 50 - 200 µL	D200	F123601	Min. 50	181,20 €
			100	
			Max. 200	
P1000 200 - 1000 µL	D1000	F123602	Min. 200	181,20 €
			500	
			Max. 1000	
P5000 1000 - 5000 µL	D5000	F123603	Min. 1000	218,40 €
			2000	
			Max. 5000	
P10 mL 1 - 10 mL	D10 mL	F161201	Min. 1 mL	218,40 €
			2 mL	
			5 mL	
			Max. 10 mL	



### Pipetman® Starter Kits

#### ◎ All-in-one savings

- Each kit contains two or three classic Pipetman pipettes plus boxes of certified-quality Gilson Diamond tips and useful accessories, such as Single™ pipette holders (except Large-volume Kit).
- Save money when compared to pipettes and accessories purchased separately.

Pipetman Starter Kits†						
Model	Part Number	Pipetman Models	Diamond Tips (racks of 96 tips)	Accessories	List price	Special price
Pipetman Micro-volume Kit	F167500	P2	Sterile filter tips	3 Single pipette holders 10 Coloris identification clips 1 Gilson Guide to Pipetting	599,00 €	†
		P10	DFL10ST and			
		P100	DF100ST			
Pipetman Starter Kit	F167300	P20	D200	3 Single pipette holders 10 Coloris identification clips 1 Gilson Guide to Pipetting	525,60 €	†
		P200	and			
		P1000	D1000			
Pipetman Large-volume Kit	F167400	P5000 P10 mL	5 x D5000 tips 2 x D10 mL tips (no racks)	2 x 10 protection filters for P5000 and P10 mL	360,00 €	299,00 €



**Special Offer!**

Purchase Pipetman Large Volume Kit for special price!\*

Offer valid October 1st – December 31st 2009

Promo Code KIT PROMO 09

Gilson maximum permissible errors are guaranteed only when Pipetman pipettes are used with the recommended Gilson Diamond® tips.

† Special offer does not apply to Pipetman Starter Kits or Microvolume Kits.